



Just Add Gas

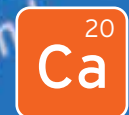

## BWB Syn Fuels Flame Photometer

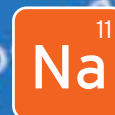

*Just Add Gas*



## Flame photometers that push the boundaries of low temperature flame analysis.

Analysis with a low temperature propane, LPG, or butane flame has been redefined by BWB Technologies.



## Our Mission

To redefine the boundaries of Flame Emission Spectroscopy (FES) by supplying unique, cost effective, reliable and accurate flame photometers (and related products) and to build strong, lasting relationships with our customers by providing the best customer service in the industry.

### About BWB Technologies

BWB Technologies is a UK company focused on the design, manufacture and sales of award winning flame photometers. Drawing on an Anglo-American team of leading industry specialists, we strive to create high quality, cost effective products that redefine what is achievable with a low temperature flame. From our manufacturing plant in Newbury, England, BWB Technologies have introduced a series of flame photometer products and accessories which exceed existing expectations in terms of specification, usability, accuracy, reporting, build quality and value for money.

### International support network

BWB Technologies are the only flame photometer manufacturer to have dedicated sales and technical support offices in Europe, The America's and the Middle East. We offer an unsurpassed level of customer support with a 24/7 support service and real people answering customer calls.

### Pushing the boundaries of flame photometry

What can be achieved with a flame photometer has been redefined by BWB Technologies. With unrivalled accuracy, low cost and ease of use the BWB flame photometer is not only the best flame photometer available but is the first real alternative to AAS (atomic absorption spectrometry) and ICP (Inductively Coupled Plasma) for measuring Lithium (Li), Sodium (Na), Potassium (K), Calcium (Ca), Caesium (Cs) and Barium (Ba). The recent introduction of industry specific products for Biological samples, Sugar, Nuclear and others further defines BWB Technologies as the leader in "next generation" flame photometers.

### Best customer service in the industry

- International telephone service available 24/7, with local access numbers in Europe, the USA and the Middle East
- An international distributor network servicing the needs of customers in over 100 countries.
- Fast delivery of spare parts and consumables backed by extensive inventory in multiple strategic locations.
- Pre-sales application specific product selection assistance
- 100's of hours of video training available in a wide choice of languages
- Technically competent customer service staff
- ISO9001 Certified manufacturer
- Excellent Globally acclaimed after-sale support and service by Value Added Distributors and Internet based assistance from GMT and PST time zones.

### Contact Information

BWB Technologies UK Ltd Unit 7, The Galloway Centre, Hambridge Lane, Newbury, RG14 5TL UK

T **Europe** + 44 (0)1787 273 451 **United States** 1 (800) 608-9870 **Middle East** +357 22 666080

E **Orders** orders@bwbtech.com **Admin** admin@bwbtech.com **Technical** technicalsupport@bwbtech.com **Sales** sales@bwbtech.com

**www.bwbtech.com** 12 months warranty

BWB Technologies UK, Ltd. warrants, subject to the conditions itemised within this document, through either BWB personnel or personnel of its authorised distributors, to repair or replace free of all charges, including labour, any part of this product which fails within the warranty time specified above, appertaining to this particular product. Such failure must have occurred because of a defect in material or workmanship and not have occurred as a result of operation of the product other than in accordance with procedures described in the instructions furnished with this product.

## Common Features

All BWB products are built around the same innovative product design and world leading technologies, sharing the same core features across the range.

### Just add gas!

BWB have revolutionised the flame photometry market by offering an instrument with all accessories supplied as standard. This means that the process of ordering a BWB flame photometer is incredibly simple. There are no complex product codes, accessory lists or ordering requirements. When you order a BWB flame photometer it will arrive with all the accessories you need.

### Up to 5 Channels of simultaneous detection and display

The ability to measure and calibrate up to 5 different ions at the same time offer huge advantages over more traditional instrumentation. As the BWB is the first and only flame photometer to offer simultaneous detection and display, our customers do not need to waste valuable time changing filters and re-calibrating when a different measurement is required.

### Built-in air compressor

Greatly reduces noise and laboratory footprint while increasing instrument stability, accuracy and reproducibility.

### Advanced calibration system

The ability to store calibration data and perform multi-mode calibrations reduces analysis time significantly and increases accuracy.

### Interfaces

Connect to your PC using USB 2. A 0-2.5volt output is supplied as standard for a chart recorder or for process control a 4-20mA output is available.

### New modern technology offers increased accuracy and reliability

Introduced in 2005, the BWB flame photometer was designed from the ground up, using new, innovative and modern technology to give unsurpassed levels of accuracy, usability and reliability. With clean lines and contemporary bright colour schemes, the instrument looks exceptionally modern and will fit into any laboratory. The instrument has a large LCD display which can be easily viewed from any angle.

### Computer connectivity as standard!

The FP-PC software package, included with every purchase, has been written with GLP compliance in mind and facilitates the creation of data and reports in .pdf and .csv formats.

### Ease of use

With exceptional ease of use, a clear backlit liquid crystal display guides users through analysis using simple menu options. The display, used in conjunction with a comprehensive membrane switch panel ensures that accurate and reliable results can be obtained with minimal training.

As modern businesses may experience a large turnover of staff, laboratory managers will appreciate the comprehensive training video supplied with every BWB flame photometer. Our training videos are also available online in over 10 languages at [www.bwbtech.com](http://www.bwbtech.com) and on BWB's YouTube Channel <http://www.youtube.com/user/BWBTechnologiesUKLtd>

### Long term savings

Excellent value for money and very low total cost of ownership.

Selective analysis Flame photometry excludes most elements which means very little interference.

### Reliable & simple

Manipulation of the flame is the only adjustment available to the operator during operation ensuring reliable and accurate measurements every time.

### Fast!

Once the instrument has stabilized it takes just a matter of seconds to obtain results.

### Collection cup

Optional for 'online', continuous measurements, the collection cup turns the BWB Flame Emission Spectrometer (FES) into an 'online', real time, 24/7 analyser with 4-20mA output linked to a single channel. Via a 2 wire output the BWB FP can seamlessly integrate with any SCADA system. Totally unique and available on all models.



# Green Credentials

We take our green responsibilities very seriously.

We constantly evaluate our processes and procedures in order to minimise the impact to the environment.

This is how we do it.

## Making our blue more green



### Our green goals

- Using LED lighting with efficient heating systems across our manufacturing facility
- Utilising only 'green' energy suppliers our electricity is supplied from 100% renewable sources and our natural gas consumption is 100% carbon neutral
- Environmentally friendly packaging created from pre-recycled content or utilising materials that can be easily recycled across the globe
- The use of paper based packaging tape to reduce our plastic consumption and allow our cardboard boxes to be easily recycled.
- Product design considerations such as recyclability and life cycle assessments are taken into account at the development phase of our new products

## Tree planting promise

We have our own in house recycling scheme for factory generated waste to ensure that our employees and the company are doing all we can to help protect and preserve the environment

In addition we have an ultimate goal to plant 1000 trees for every company employee by 2025

We're proud to be the first Flame Photometer manufacturer to plant trees for every instrument sold. We're committing to planting at least 50 trees for every model of Flame Photometer we sell from 01 May 2021



Striving to meet as many of the United Nations 17 sustainability development goals is a key priority for us in our commitment to protect the environment. We have adapted our processes and products to meet the following goals:



# BWB SYN FUELS

## Reliable and accurate analysis of all four atoms of interest

This powerful and growing global industry now has the attending facilities for Synthetic Fuels or Bio-Fuels production. For this derivative product BWB Technologies offers a hybridized Flame Photometer we have named Syn Fuels.

The BWB SYN FUELS Flame Photometer is the first and only 4 channel flame photometer with simultaneous detection and display of all 4 elements. A high quality and high performance instrument, the BWB SYN FUELS employs modern technology to improve both accuracy and stability while significantly reducing analysis time.

The fluid handling sections of the BWB SYN FUELS Flame Photometer incorporate proprietary engineered thermoplastics that are highly chemically resistant to withstand the solvents encountered in Synthetic Fuels production facilities.

The majority of our Global clients purchase a BWB Sugar FES for the beet or cane refining plant and the Syn Fuels FES for the bio-conversion plant site.

For non Sugar processors embracing Syn Fuels, the BWB XP is typically the BWB FES of interest during the bio-conversion processes.



### + FEATURES

- SIMULTANEOUS DETECTION AND DISPLAY OF ALL 4 ATOMS OF INTEREST
- "IRS" (INTERNAL REFERENCE STANDARD) AVAILABLE
- INTUITIVE USER INTERFACE FOR TRUE EASE OF USE
- DISPLAY PROMPTS STEP BY STEP OPERATION
- BUILT-IN AIR COMPRESSOR
- SOLUTIONS AND LABWARE INCLUDED
- DATA SHARING VIA PC LINK
- OPERATOR INDEPENDENT DETERMINATION OF RESULTS
- 4 USER SELECTABLE UNITS OF MEASURE
- USER SELECTABLE DECIMAL PLACES
- OPTIONAL INTEGRATED PRINTER USES READILY AVAILABLE PAPER
- IQ, OQ, PQ WEB-BASED CERTIFICATION AVAILABLE
- CORRECTION OF CA FOR THE INTERFERENCE FROM HIGH LEVELS OF NA



*Just Add Gas*

# Technical Data

*Just Add Gas*

### Sample rate

2 - 3.5ml/min

### Tubing materials

PTFE and Teflon

### Required desk space for instrument

50cm (H) x 45cm (W) x 45cm (D)

For Safety reasons the Flame photometer requires 1m of unobstructed space above to allow dissipation of heat from the chimney.

### Instrument size

Syn Fuels Flame Photometer

51cm (H) x 39cm (W) x 35cm (D)  
(20in x 15in x 14in)

### Shipping

65cm (H) x 49cm (W) x 43cm (D)  
(25in x 19in x 17in)

### Weight

Syn Fuels Flame Photometer - 14kg (30.8lbs)

Shipping - 25kg (55.1 lbs)

### Optimal range

Single point calibration

Na 0.1 - 60ppm

K 0.05 - 100ppm

Li 0.05 - 50ppm

Ca 1.0 - 100ppm

### Multi Point Calibration

Na 0.1 - 1000ppm

K 0.05 - 1000ppm

Li 0.05 - 1000ppm

Ca 1.0 - 1000ppm

### Reproducibility

<1% Coefficient of variability for 20 consecutive samples over ten minutes at concentrations of 100ppm or less. (After instrument stabilisation).

### Specificity

Na/K/Li = <1% to each other when equal in concentration at <100ppm

### Limit of detection (LOD) and limit of quantification (LOQ) displayed

LOD	LOQ
Na - 0.03ppm	Na - 0.1ppm
K - 0.02ppm	K - 0.05ppm
Li - 0.02ppm	Li - 0.05ppm
Ca - 0.3ppm	Ca - 1ppm

### Time to stability

Less than 15 seconds after sample is introduced into the flame

### Ions measured

Simultaneous measurement of Na, K, Li, Ca

### Interfaces

USB

0-2.5 volt output (based on sample concentration linked to element of users choice)

Optional 4-20mA output in place of the above

Optional integrated printer

.CSV and .PDF generated reports and files via FP-PC software

Collection Cup fitted as standard allowing live reading

### Recommended minimum warm up time based on ambient temperature

21°C - 25 minutes

### Power requirements

100V - 250V AC at 50 or 60Hz

### Fuel requirement

Propane, Butane or Natural Gas\* regulated to 19Bar. Flow rate of 0.4l/min

\*with modifications. BWB Technologies recommends either Propane or Butane for optimum results.

### Readout

LCD, four line, alpha numeric, back lit.

## + WHAT'S IN THE BOX

- Power lead to suit region
- Gas hose
- Aspiration kit
- Manuals covering all aspects of the Syn Fuels FP
- Quick start guides
- Warranty registration form
- Certificates of analysis for all provided fluids
- Material safety data sheets for all provided fluids
- The FP-PC software installation USB
- 150ml of BWB recommended cleaning solution
- 150ml calibration fluids for all 4 elements at 10,000ppm.
- 100ml conical flask with screw lid
- 100 disposable sample cups
- 1L of de-ionised water
- 2m of waste tube
- 10 x 10ml pipettes
- 10 x 1ml pipettes



Proudly Engineered  
and Manufactured  
in England

[www.bwbtech.com](http://www.bwbtech.com)

Manufactured in Newbury, County of Berkshire, England.  
Registered office: BWB Technologies Ltd Unit 7, The Galloway Centre, Hambridge Lane, Newbury, RG14 5TL UK

## FP-PC and BWB Drive

### FP-PC app

The BWB flame photometer has been configured to easily satisfy the most sophisticated IT requirements. The FP-PC app provides a useful interface to the BWB Flame Photometer (FP) facilitating the following operations:

- Normal calibration Report production
- Logging of streaming data for continuous on-line applications. (A 4-20 mA signal can also be output proportional to a calibrated channel on the FP).
- Extraction of data from the flame photometer when operated in automatic or semi-automatic mode with or without auto-sampler / auto-diluter known as Automatic Fluid Handling System (AFHS)

Calibration manipulation including:

1. Upload from / download to the FP
2. Storage / retrieval of calibrations using the PC or system hard drive

These functions are provided by the Welcome, Calibrate, Report, Logging, AFHS tabs together with one more (XP-Cal or BIO-Cal) that is dependent on the mode that the FP is being operated in.

The FP-PC app, included with every purchase, has been written with GLP compliance in mind and facilitates the creation of data and reports in PDF format.

### BWB drive

Continuing to building on modern technologies, BWB have created an industry leading USB flash drive packed full of features and resources designed to aide your daily analysis and BWB experience. Exclusively for BWB customers, this self-contained resource includes:

- Tutorial videos in 10 languages
- PC and Mac compatibility
- FP-PC app installation package (PC only)
- Links to BWB's growing online video library (internet connection required)
- Manuals, guides, warranty information and registration PDF
- A list of international support and distributor contact details
- 13 GB of free space for you to use as you please!

This wealth of information is provided in a beautifully designed BWB branded USB memory stick, which includes a built in flash light!

### [www.bwbtech.com](http://www.bwbtech.com)

The constantly evolving [www.bwbtech.com](http://www.bwbtech.com) maintains the companies striking bold colours and clean lines ensuring a crisp and easy to use interface packed full of new and exciting resources. With weekly updates the BWB website is at the forefront of the industry when it comes to education and resourcefulness, among many features it includes:

- FAQs
- FP Practical Applications
- User guides and reference materials
- Product specs and technical documents
- A vast video support library
- Support documents
- Local distributor lists
- Tech Support and Sales
- Support contact information



### + THANK YOU

BWB takes pride in working closely with its distributors from across the globe. Spearheaded by our New Media Consultant Simon Hustings, The BWB flash drive would not have been possible without the direct contribution from some of our amazing distributors.

Thank you to: Wolfgang, Mariusz, Andre, Jereon, Nhung, Stonepepple, Kaoriko, Fernando, Alvaro, Alexandra, Dr. Muhamed, Dimitri, Francois and Thiemdej.



## **IQ OQ PQ Certification**

### **IQ OQ PQ is a unique validation and certification program from BWB.**

The programme has been specially created to ensure the end user has the technical knowledge required to install and initialise their BWB Flame Photometer and the technical know-how and confidence to operate any of our 2015 and onwards range of flame photometers correctly and efficiently on a daily basis.

#### **IQ - Installation Qualification**

The installation qualification covers in detail the key steps required to correctly set up the BWB range of photometers and Auto Fluid Handling System. IQ stages include: maintenance and service plans, the installation environment and health and safety requirements.

#### **OQ - Operation Qualification**

The operation qualification covers the steps required to ensure that the BWB flame photometer and AFHS are operating correctly.

The user will undertake a series of hardware checks and test readings, providing the results to BWB for analysis.

#### **PQ - Performance Qualification**

The performance qualification has been designed to ensure that the instruments perform to the original factory specifications.

Once the IQ and OQ sections have been completed, the user runs a sample manually or using AFHS with the solutions, prints out the results and sends to BWB for analysis.

After completion of the qualification and confirmation from BWB, a signed certificate will be issued to the user.

IQ OQ PQ is a unique and effective program built to ensure technical ability, confidence and peace of mind for all users of the BWB 2015 and onwards Series of Flame Photometers and Auto Fluid Handling Systems.

To find more about IQ OQ PQ get in touch at [info@bwbtech.com](mailto:info@bwbtech.com)

### **+ THANK YOU**

Over the years, BWB Technologies has forged strong ties within the industry and IQ OQ PQ certification is the evolution of those bonds. A big thank you to Wolfgang, Alexandra, "Eddie", Fernando, Iveta and Anna, Team Leader Harry R. Libengood - Company Co-Founder, and Hozan Edwards for their focus and dedication in developing IQ OQ PQ.

## Testimonials



**Alvaro Apéstegui, Ph.D., Laboratorios Labiclin, Costa Rica** **BWB BIO**

"I have been working with ion selective instruments for a while, but after numerous errors with the lipemic samples, and the high running and maintenance costs associated with the instrument I started looking for a new flame photometer. I needed an instrument that gave me good precision, repeatability and an instrument that is viewed as the gold standard for ions such as sodium, potassium and lithium.

I have analysed many instruments but they always had something that was not adequate for my needs. I decided to go for a BWB BIO because it had very good construction and offered up to five different ions with calibration curves of up to fifteen points.

I have had the instrument for almost a year and must say that it is here to stay because of its sturdy construction, really fantastic electronics and most importantly, it will not allow you to make mistakes. Being able to actually see the stability of the readings and having calibrations that are very stable is comforting.

The communication with the people from BWB Technologies is fluent and professional. I have never had a question unanswered and responses normally come in less than 12 hours which says a lot since I am in Costa Rica and they are in England!

I strongly recommend this instrument, and now that it can have automated sampling, makes it even easier and equals the automatic capabilities of the ion selective instruments."



**Gregory Bathfield, Electrical & Control Specialists Ltd, Mauritius** **BWB SUGAR**

"For the sugar industry, the flame photometer is used to measure Potassium which is directly related to sugar. The flame photometer is thus used to detect sugar in condensate return to the boiler. The BWB photometer is more user friendly and easier to use than others we came across. It is really good value for money."



**Wolfgang Glock, THG-Glock Instruments, Germany** **BWB XP**

"My experience I have made with the BWB-XP Flame Photometer is, that it is easy to install, easy to handle, gives accurate measurements and with the software the operator has a better summary of the work he has done in the Lab. The excellent service of BWB- Technologies and the excellent performance of the instrument makes the BWB- XP Flame Photometer the "best Flame Photometer which is on the market". At the moment no competitor can beat this instrument in performance and price."



**Richard Griffiths, Fujifilm Imaging Colorants Ltd, England** **BWB XP**

"We are very happy so far with the purchase of our BWB-XP flame photometer. We are using the unit for a colours application and had upgraded from a older Sherwood model. I carry out 5 point calibrations on each ion and these calibrations appear to hold from first thing in the morning until the end of the day. The BWB-XP produces accurate values and the ability to measure 5 elements simultaneously is less time consuming than the older instrument and has, on occasion, flagged up the presence of ions which we were not expecting to see and may not have been looking for on the older instrument."

**DELIGHTED WITH YOUR BWB FLAME EMISSION SPECTROMETER? EMAIL US YOUR TESTIMONIAL AND WE WILL HAPPILY RESPOND WITH 5 BWB USB FLASH DRIVES! A GREAT KEY RING HOLDER, STORAGE SPACE AND TORCH FOR THE ENTIRE FAMILY!**

# Distributors

We have a global network of first class distributors.  
Visit our website, [www.bwbtech.com](http://www.bwbtech.com) for full details.





BWB Technologies Ltd Unit 7, The Galloway Centre, Hambridge Lane, Newbury, RG14 5TL UK

T **Europe** + 44 (0)1787 273 451 **United States** 1 (800) 608-9870 **Middle East** +357 22 666080

E **Orders** orders@bwbtech.com **Admin** admin@bwbtech.com **Technical** technicalsupport@bwbtech.com **Sales** sales@bwbtech.com

**www.bwbtech.com**

BWB Technologies Ltd. warrants, subject to the conditions itemised within this document, through either BWB personnel or personnel of its authorised distributors, to repair or replace free of all charges, including labour, any part of this product which fails within the warranty time specified above, appertaining to this particular product. Such failure must have occurred because of a defect in material or workmanship and not have occurred as a result of operation of the product other than in accordance with procedures described in the instructions furnished with this product.

Manufactured in Newbury, County of Berkshire, England.  
Registered office: BWB Technologies Ltd, Unit 7 Galloway Centre, Newbury, Berkshire, RG14 5TL, UK.



**Proudly Engineered  
and Manufactured  
in England**